

# iUL

# iPeak®+

## Lateral Flow Reader

### Main Benefits

#### Connect & integrate!

Transmit data to portable devices such as mobile phones or tablets, cloud, LIS, LIMS, HIS, EMR. Connect the reader to local networks, or the Internet, iPeak®+ meets the requirements of different environments of use, guaranteeing the security of your data. Know more about iPeak®'s connectivity here.

#### Utmost sensitivity

Detect and even quantify your target molecules even if they are in low concentrations thanks to the Flash Eye Technology.

#### Track the sample, track the results!

Equipped with an embedded barcode reader, read barcodes, DataMatrix, or any code on virtually any medium - paper, plastic cards, mobile phones, or LCD displays.

#### Stand alone

Portable and lightweight, iPeak®+ is ideal for any kind of environment.

#### User-friendly

4.3" capacitive and intuitive screen. It is essentially useful for non-trained professionals.

#### Compliance

CE marked. IUL sets out the criteria for a quality management system certified ISO 9001 and ISO 13485:2016 Medical devices - Quality management systems.

# MONLAB

SERVICES FOR CLINIC DIAGNOSTIC LABORATORIES

### Product Overview

iPeak® expands its family with the iPeak®+, colorimetric lateral flow reader for your large cassettes, up to 7 lines in one strip and up to 5 strips in one cassette.

With the breakthrough Flash Eye Technology, iPeak®+ is the perfect match for your lateral flow tests, enhancing the colors of the strip, and removing the subjectivity of the human eye.

### When optics matter


Point-of-care | Health & Veterinary Diagnostics | Food Safety | Drug testing | Environmental | Research



## Flash Eye Technology

Flash Eye Technology is the combination of excellence in image acquisition and processing and the IUL's know-how. Flash Eye technology enhances the colors of any colorimetric lateral flow test with an outstanding dynamic range far qualitative, and quantitative assays.

iPeak®+ brings a camera that captures the image of the test illuminated from LED lights situated in the most studied geometry, delivering precise and uniform lighting. The Flash Eye Technology allows to capture the best image of the test making the sensitivity of the reader its cornerstone. Know more about this technology here.

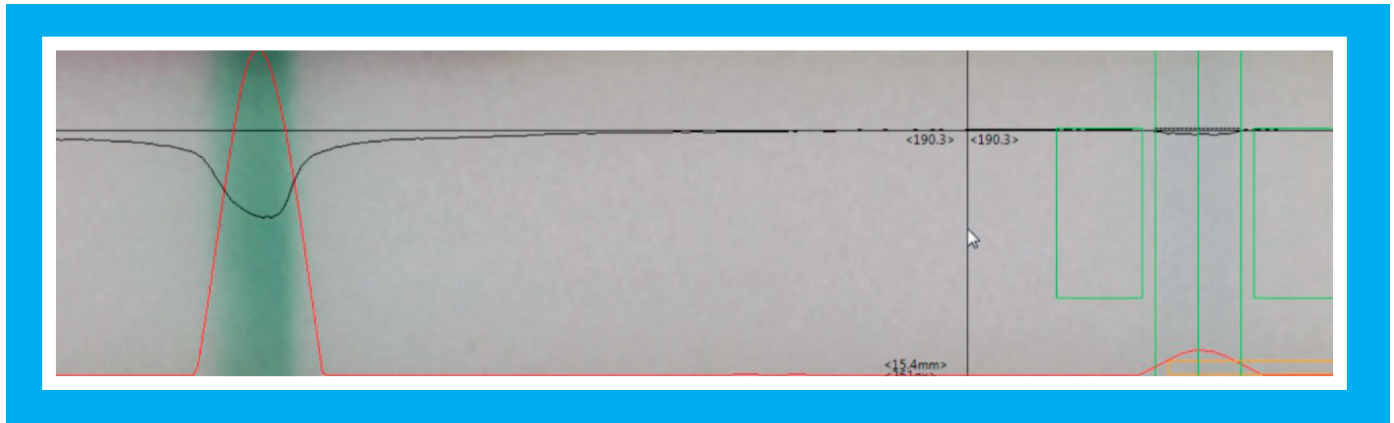

## Rapid tests available

An affordable and efficient system for an automated lateral flow reading.

The iPeak® Lateral Flow reader can help you certify your system for rapid test reading.

## Nasopharyngeal tests

Respiratory Adenovirus

RSV

Influenza A+B

Strep A

Covid19 AG

## Faeces TESTS

C. difficile A/B

C. difficile GDH

C. difficile GDH + A/B

Rota Adenovirus

Giardia

Crypto-Giardia-Entamoeba

Campylobacter

Norovirus

H.pylori

## Urine TEST

Legionella

S. pneumoniae

Legio-Pneumo Combo

Hcg

S. pneumoniae (strips)

Individual drugs of abuse

Multiple drugs of abuse (iPeak+)

## Serum /Whole blood tests

Covid19 Ac IgG/IgM



## Performance and Specifications

### Performance

Internal memory	16 Gb
Stored results	Results, raw, and processed images. Depending on the internal memory
Number of test lines per strip	Up to 7 lines / strip Up to 5 strips
Measurement time	15 seconds, one strip, 2 lines
Technology	Colorimetric lateral flow reader with the IUL's Flash Eye Technology camera - based
Languages	Español, English, Deutsch, Français, Italiano, Polski, تخبر علما

### Technical Specifications

Connectivity	16 Gb
Security & Exportation	HL7 v2.8 compliant. Connections https with SSL TLSv1 TLSv1.1 TLSv1.2 TLSv1.3. Authentication: transport level 2 with token key management. Also can export to an external SFTP server or directly from a slim USB
Audit Trail	HL7 v2.8 compliant. Connections https with SSL TLSv1 TLSv1.1 TLSv1.2 TLSv1.3. Authentication: transport level 2 with token key management
Sample traceability	Datamatrix, QR, and barcode
Battery	4 h autonomy
Power range	100-240 V
Frequency range	50-60 Hz
Maximum power	18 W
Air temperature	15-35°C
Relative humidity	15-75%
Place of operation	For indoor use only
Maximum dimensions of the cassette (LxWxH)	Up to 80 x 56 x 6 mm / 3.1 x 2.2 x 0.2 in
Dimensions (LxWxH)	150 x 140 x 185 mm / 5.9 x 5.5 x 7.2 in
Weight	900 g / 1.98 lb
Screen	Color, capacitive, 4.3" touchscre



### Order information

- iPeak®+ CAT.  
No. 900034000
- Self-adjusting light drawer  
CAT. No. 900012008
- Ticket printer (Epson)  
CAT. No. 90003069
- Cable for the printer (Epson)  
CAT. No. 90004874



### Associate Distributor

Domaine de Pissaloup  
16, rue d'Alembert  
Buroplus HTA2  
78190 TRAPPES

T. : +33(0)1 30 68 80 00

F. : +33(0)1 30 68 80 02



#### MONLAB S.L.

Selva de Mar, 48  
08019 Barcelona - Spain

Tel. +34 93 433 58 60

mn.sales@monlab.es - www.monlab.com

### Video

Scan the QR code  
and watch the video

